

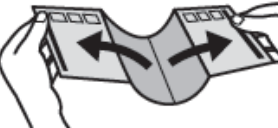











LGG Study Stool Sample Instructions

- Do not open the kit until after you are in the restroom.
- Do not collect your sample until you are sure you can drop the sample in the mailbox within 24 hours.

Procedure:	
	<p>1. Two toilet accessories are provided in the event that you are unable to collect a sample on the first attempt. Caution: Do NOT allow toilet water, urine, detergent or fragrance to come in contact with the toilet accessory.</p>
	<p>2. Empty your bladder before beginning the collection. Collect fecal sample free of urine or toilet water. Toilet paper or tissues may be required.</p>
	<p>3. Carefully peel open edge with ▲; repeat for edge with ■.</p>
	<p>4. Attach adhesive surface of accessory to BACK of toilet seat with adhesive close to the outer edge of the seat and press firmly.</p>
	<p>5. While holding the yellow tube top, unscrew ONLY the purple cap and set aside for later use.</p> <p>Important: Do NOT remove the yellow tube top. Do NOT spill the stabilizing liquid in the tube.</p> <div style="text-align: right;">  </div>
	<p>6. Use the spatula to collect a small amount of fecal sample.</p> <div style="text-align: center;">  <p>Actual size of fecal sample.</p> </div>
	<p>7. Transfer the fecal sample into the top portion of the yellow tube top. Repeat until the sample fills the yellow tube top.</p> <p>Important: Do NOT push sample into the tube.</p> <div style="text-align: right;">  </div>
	<p>8. Scrape horizontally across the tube top to level the sample and remove any excess.</p> <p>Wipe exterior of tube and top with toilet paper or tissue as needed.</p>

Procedure (continued):	
	<p>9. Pick up the purple cap and screw onto the yellow tube top until tightly closed.</p> 
	<p>10. For a minimum of 30 seconds, shake the sealed tube as hard and fast as possible in a back and forth motion.</p>
	<p>11. The fecal sample will be mixed with the stabilizing liquid in the tube; not all particles will dissolve.</p> <p>Important: Continue shaking if large particles remain as shown in Figure A.</p>
	<p>12. Drop used accessory into toilet. Wait 1 minute for paper to become soft, then flush. Alternatively, discard in garbage.</p>
	<p>13. Place spatula in original packaging or wrap in toilet paper and discard in garbage.</p> <p>Important: Send the sample for processing following the delivery instructions supplied separately by the kit provider.</p>

After you have collected your sample:

- Record the date and time you collected your sample on the outside of the tube.
- Place the sample in the self-addressed envelope provided to you.
- Seal the envelope.
- Drop the envelope off at the nearest mailbox or post office.